

Introducing the new

Davies® Gold Series dermatome and electric motor system

The Davies® Gold Series range of dermatomes are precision engineered instruments with a proven track record for quality and dependability. Two versions are available, the Simplex with a 48mm cutting width and Duplex with a 73mm cutting width. Both dermatomes are designed for use with the CM2 controller and Micro motor, however they can be driven by alternative motor systems. - Please ask for further details.

- Compact and lightweight
- Wide range of cutting depths
- Selection of cutting widths
- Robust, proven design
- Dermatome heads, motor unit and cable are fully sterilisable
- Very accurate and consistent cutting performance
- Powerful, variable speed, high torque motor unit
- CE mark conformance
- Easily maintained and serviced



Please contact De Soutter Medical Limited for further details of the Davies® Gold Series range of dermatomes

Dermatome Specification

Model	Simplex	Duplex
Cutting width (mm)	48	73
Cutting thickness range-7 settings	0.08 - 0.38 mm	
Operating speed range (RPM)	2000 - 15000	
Weight (Kgs)	0.25	0.32

Micro motor and power supply

Both dermatomes quickly fit onto the universal micro motor which is powered by the variable speed CM2 controller. The motor can be controlled either by the footswitch or directly from the controller itself. Please specify voltage and plug configuration when ordering.

Power supply specification

Input voltage:	230 (120) VAC/50 (60) Hz
Power:	130VA
Dimensions (mm):	100 x 85 x 165
Weight (Kgs):	1.68/0.12 controller/motor


Optional

Reduction sole plates	Part No.
Simplex to 23 mm cutting width	8530
Duplex to 48 mm cutting width	8540
Duplex to 23 mm cutting width	8550



Sterilisation container

This purpose designed compatible delivery system will ensure your instruments are maintained in good condition at all times. The container is 100% sterilisable by steam, dry heat ('flash') to 160°C, ethylene oxide (ETO) or gamma irradiation.

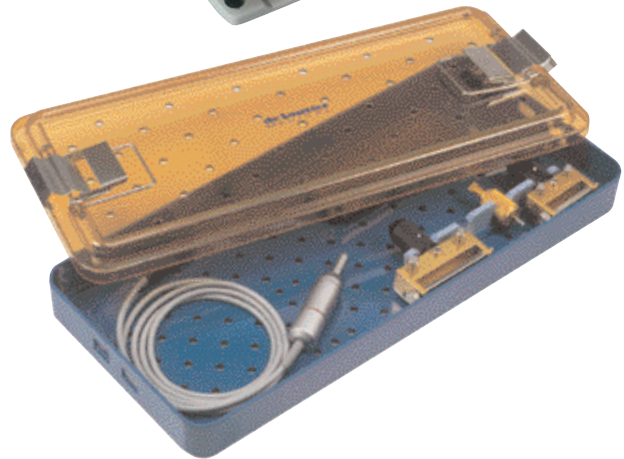

How to order

Simplex system - 48 mm cutting width

System components	230V - UK plug	230V - Euro plug	120V - US plug
GS-102 Simplex head complete	222884	222884	222884
CM2 motor controller and footswitch	1280884	1280994	1281004
Simplex blades (pack of 5)	S11-102	S11-102	S11-102
Sterilisation container	8570	8570	8570
Complete system ordering code	GSA-102/230UK	GSA-102/230EU	GSA-102/120US

Duplex system - 73 mm cutting width

System components	230V - UK plug	230V - Euro plug	120V - US plug
GD-103 Duplex head complete	222994	222994	222994
CM2 motor controller and footswitch	1280884	1280994	1281004
Duplex blades (pack of 5)	S11-103	S11-103	S11-103
Sterilisation container	8570	8570	8570
Complete system ordering code	GDA-103/230UK	GDA-103/230EU	GDA-103/120US

de Soutter
MEDICAL

De Soutter Medical Limited

River Park, Billet Lane, Berkhamsted
Hertfordshire HP4 1HL, U.K.

Tel: 01442 860300 Fax: 01442 860333

Email: info@de-souttermedical.com

Web: http://www.de-soutter.com

